



Autoclaves for sterilization "Presoclave III" 50, 80 and 150 litres.

MICROPROCESSOR TEMPERATURE AND CYCLE CONTROL AND ATMOSPHERIC PURGE.
ELECTRONIC CONTROL OF TIME AND TEMPERATURE RANGE FROM 115 °C TO 134 °C. (0.62 TO 2 bar).
CONFORMS TO APPARATUS PRESSURE REGULATION DIRECTIVES.
VERTICAL MODEL: 50, 80 AND 150 LITRES.



SAFETY:


- SAFETY OVER PRESSURE VALVE:** PREVENTS OVER PRESSURE ABOVE THE MAXIMUM SET LIMIT.
- SAFETY PRESSURE SWITCH:** DISCONNECTS POWER TO THE HEATER IF THE UNIT EXCEEDS THE MAXIMUM PRESSURE LIMIT.
- SAFETY THERMOSTAT AND WATER LEVEL:** DISCONNECT POWER TO THE HEATER IF INSUFFICIENT WATER OR OVER TEMPERATURE.
- SAFETY LID:** PREVENTS THE LID BEING OPENED IF THE CHAMBER IS UNDER PRESSURE.
- SAFETY INSULATION:** THERMALLY INSULATED LID .

An economic alternative, semi-automatic

FEATURES

- External case made from stainless steel.
- INCOLOY stainless steel heater elements, hard waring, resistant to corrosion and high temperatures.
- Sterilization chamber and door made from stainless steel.
- Automatic purge to atmosphere
- Lid, 6 point closure.

CONTROL PANEL

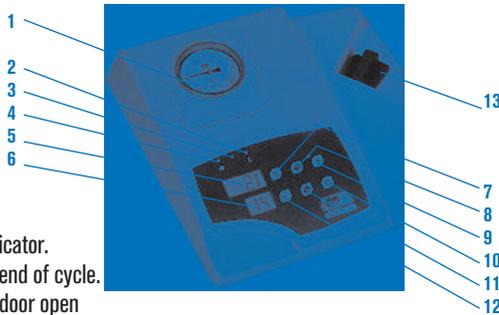

1. Pressure gauge.
2. Illuminated alarm indicator.
3. Illuminated indicator end of cycle.
4. Illuminated indicator door open
5. Temperature display.
6. Time display
7. Temperature set pushbutton (range from 115 °C to 134 °C).
8. Pushbutton increase value.
9. Push button START
10. Push Button STOP
11. Push button Reduce value.
12. Push button select time, (range from 3' to 99').
13. Push button open lid .



MODELS	Part No.	Capacity litres	Ø	Height (usables cm)	Height / Width / Depth (exterior cm)	Power W	Weight Kg
PRESOCLAVE III 50	4001758	50	30	70	117 48 62	3600	75
PRESOCLAVE III 80	4001759	80	40	60	107 58 72	3600	100
PRESOCLAVE II 150	*4002137	150	50	70	118 80 95	7500	225

Does not include baskets and drums, order separately. Part No. 4002137 due to its consumption is manufactured for three-phase current at 230 V or 400 V to be specified.

ACCESSORIES: AISI 304 Stainless steel wire baskets and drums.


- For Presoclave III 50, Part. no: 4001758.
Capacity, 3 baskets Ø 25 and 20 cm. height. Part No. **1000495**
- For Presoclave III 80, Part. no: 4002136.
Capacity, 2 baskets Ø 36 and 28 cm. height. Part No. **1000496**
- For Presoclave II 150, Part. no: 4002137.
Capacity, 2 baskets Ø 44 and 22 cm height. Part No. **1000780**
- Sterilization drum made from AISI 304 stainless steel.
For Presoclave III 50, Part. no: 4001758.
Capacity, 3 drum Ø 25 and 18 cm. height. Part No. **1002418**
- For Presoclave III 80, Part. no: 4002136.
Capacity, 2 drum Ø 35 and 27 cm. height. Part No. **1002419**

Detergent Bio-Sel.

Water distillation unit especially for autoclaves.

USB adapter model. Pen-Drive included (Memory board) for data storage.

Part No: **4120131**

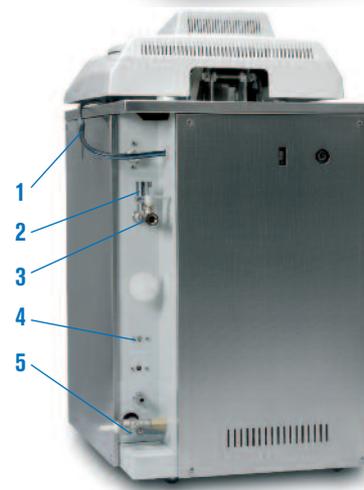

BACK VIEW:

1. Mains power.
2. Safety valve.
3. Steam empty valve.
4. Safety thermostat
5. Drain selector valve.

COMPLEMENTS That needs to be factory installed:

- Pt100 second sensor only temperature reader. Part No. **4001219** (see page 63).
- Printer of continuous paper roll with print out of temperature, time and mode. Part No. **4120113** (see page 64).
- Loading and unloading crane. Part No. **1002426** (see page 70).